



QuickVue®  
Dipstick Strep A TEST

**CLIA Complexity: WAIVED**

For *in vitro* diagnostic use.

**Rx ONLY**

A symbols glossary can be found at [quidel.com/glossary](http://quidel.com/glossary).



### INTENDED USE

The QuickVue Dipstick Strep A Test is intended for the rapid, qualitative detection of Group A Streptococcal antigen from throat swabs or confirmation of presumptive Group A Streptococcal colonies recovered from culture. The test is to be used to aid in the diagnosis of Group A Streptococcal infection. For use by healthcare professionals.

### SUMMARY AND EXPLANATION

Group A Streptococcus is one of the most important causes of acute upper respiratory tract infection. Early diagnosis and treatment of Group A Streptococcal pharyngitis has been shown to reduce the severity of symptoms and serious complications such as rheumatic fever and glomerulonephritis.<sup>1</sup> Conventional identification procedures for Group A Streptococcus from throat swabs involve the isolation and subsequent identification of viable pathogen techniques that require 24 to 48 hours or longer for results.<sup>2</sup>

### PRINCIPLE OF THE TEST

The QuickVue Dipstick Strep A Test is a lateral-flow immunoassay utilizing Quidel's patented antibody-labeled particles. The test detects either viable or nonviable organisms directly from throat swabs or culture colonies within 5 minutes.

To perform the test, a Throat Swab specimen is collected. Antigen is extracted from the Swab specimen with Reagents A and B. The Dipstick is then added to the extracted sample.

If the sample contains Strep A antigen, a pink-to-red Test Line along with a blue procedural Control Line will appear on the Dipstick, indicating a positive result. If Strep A antigen is not present, or present at very low levels, only a blue procedural Control Line will appear.

### REAGENTS AND MATERIALS SUPPLIED

***Each kit contains:***

- Individually packaged Dipsticks (25 or 50): Dipsticks coated with rabbit polyclonal anti-Group A Streptococcus
- Extraction Reagent A (1): Contains 4 M sodium nitrite
- Extraction Reagent B (1): Contains 0.2 M acetic acid
- Sterile Throat Swabs (25 or 50)
- Tubes (25 or 50)
- Positive Control (1): Heat-inactivated Group A Streptococcus with 0.2% sodium azide

- Negative Control (1): Heat-inactivated Group C Streptococcus with 0.2% sodium azide
- Package Insert (1)
- Procedure Card (1)

## WARNINGS AND PRECAUTIONS

- For *in vitro* diagnostic use
- Do not use beyond the expiration date printed on the outside of the box.
- The Dipstick must remain sealed in the protective foil pouch until just prior to use.
- Reagent B contains an acidic solution. If the solution contacts the skin or eye, flush with large volumes of water.
- Use of nitrile or latex gloves is recommended when handling the Extraction Reagents within this kit.<sup>3,4</sup>
- Do not interchange reagent bottle caps.
- If Reagent B is green prior to mixing with Reagent A in the Tube, do not use and contact Technical Support.
- To obtain accurate results, you must follow the Package Insert instructions.
- Testing should be performed in an area with adequate ventilation.
- Dispose of containers and unused contents in accordance with Federal, State and Local regulatory requirements.
- Wear suitable protective clothing, gloves, and eye/face protection when handling the contents of this kit.
- Wash hands thoroughly after handling.
- For additional information on hazard symbols, safety, handling and disposal of the components within this kit, please refer to the Safety Data Sheet (SDS) located at [quidel.com](http://quidel.com).

## KIT STORAGE AND STABILITY

Store the kit at room temperature, 59°F to 86°F (15°C to 30°C), out of direct sunlight. Kit contents are stable until the expiration date printed on the outer box. Do not freeze.

## SPECIMEN COLLECTION AND STORAGE

Collect throat swab specimens by standard clinical methods. Consult standard reference procedures such as the collection method described by Miller and Holmes.<sup>5</sup> Depress the tongue with a tongue blade or spoon. When swabbing the throat, be careful not to touch the tongue, sides or top of the mouth with the Swab. Rub the Swab on the back of the throat, on the tonsils and in any other area where there is redness, inflammation or pus.

Use rayon tipped swabs to collect throat specimens. Do not use calcium alginate, cotton tipped or wooden shaft swabs.

It is recommended that swab specimens be processed as soon as possible after collection. Swabs can be held in any clean, dry plastic tube or sleeve up to 72 hours at room temperature (15°C to 30°C), or refrigerated (2°C to 8°C) before processing. The use of charcoal or agar medium is not recommended.

If a culture is desired, lightly streak the swab on a 5% sheep blood agar plate before using the Swab in the QuickVue Dipstick Strep A Test. Do not perform the QuickVue Dipstick Strep A Test before streaking the Swab, as the Extraction Solution will destroy the bacteria on the swab, thereby rendering the organism incapable of successful culturing. Alternatively, throat swab specimens can be obtained by dual Swabs or by two sequential Swabs for the culture procedure.

## CULTURE CONFIRMATION

The QuickVue test can be used to confirm the identification of Group A Streptococcus on blood agar plates. Lightly touch a colony using a sterile Swab. Do not sweep the plate. Follow the instructions in the TEST PROCEDURE section to test the Swabs.

## QUALITY CONTROL

### *Built-in Control Features*

The QuickVue Dipstick Strep A Test provides three levels of internal procedural controls with each test run. For daily quality control, Quidel recommends documenting that these internal controls were checked for the first sample tested each day.

- The color of the Extraction Reagent changes from clear to green as the reagents are mixed together. The color change is an internal extraction reagent control and is an indication that the reagents are mixed and functioning properly.
- The appearance of a blue Control Line is an internal control. The Dipstick must absorb the proper amount of sample and the Dipstick must be working properly for the blue Control Line to appear. Additionally, the appearance of the Control Line indicates that capillary flow occurred.
- A clear background is an internal background negative control. If no interfering substances are in the sample and the Dipstick is working properly, the background in the Result area should be white to light pink within 5 minutes and not interfere with the reading of the test result.

### *External Quality Control Testing*

External controls may also be used to demonstrate that the reagents and assay procedure perform properly.

Quidel recommends that Positive and Negative Controls be run once for each untrained operator, once for each new shipment of kits - provided that each different lot received in the shipment is tested - and as deemed additionally necessary by your internal quality control procedures, and in accordance with Local, State, and Federal regulations or accreditation requirements.

If the controls do not perform as expected, repeat the test or contact Quidel Technical Support before testing patient specimens.

## ASSAY PROCEDURE


- Do not remove Dipsticks from the foil pouch until ready to perform the assay.
- To avoid cross contamination, do not allow the tip of the reagent bottles to come in contact with sample Swabs.

## TEST PROCEDURE

### ***Important:***

- Gloves should be worn when handling samples.
- Do not use Reagent B if the solution is green prior to mixing with Reagent A in the Tube. If this occurs, contact Technical Support.


1. Just before testing, add **3 DROPS** of Reagent A and **3 DROPS** of Reagent B into a clean Tube. This solution should turn **green**.

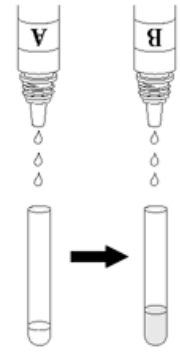


When adding drops, hold bottle vertically so that a complete drop forms.

2. Immediately add the patient Swab sample to the Tube. **Squeeze** the bottom of the Tube so the Swab head is compressed. Rotate the Swab a **minimum of 5 times**.

**Keep swab in Tube for 1 minute.**

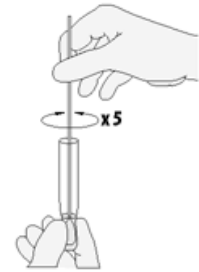

3. Express **all** liquid from the Swab against the inside of the tube. **Squeeze** the Swab firmly as it is removed from the Tube. Discard the Swab.

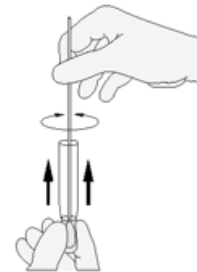

4. Remove the Dipstick from the foil pouch. Place the Dipstick into the Tube with the arrows of the Dipstick pointing down. Do not handle or move the Dipstick until the test is complete and ready for reading.



5. Read result at 5 minutes. Some positive results may appear sooner.



---

## INTERPRETATION OF RESULTS

### Positive Result:


Any pink to red Test Line along with any shade of a blue procedural Control Line is a positive result for the detection of Group A Streptococcus antigen.

### Negative Result:

A blue procedural Control Line and no pink Test Line is a presumptive negative result.

### Invalid Result:

The test result is invalid if a blue Control Line is not visible at 5 minutes. If this occurs, retest using a new sample and a new Dipstick or contact Technical Support.

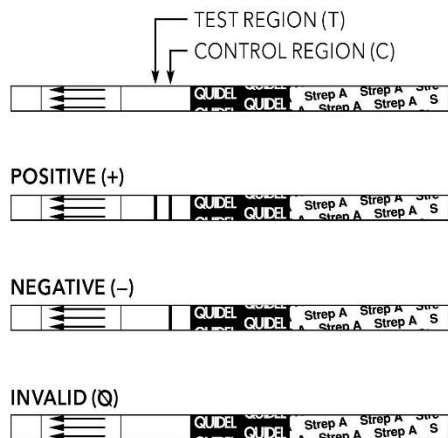

## QC TESTING PROCEDURE

- Follow the instruction procedures in the TEST PROCEDURE to dispense the Extraction Reagents into the Tube (step 1).
- Vigorously mix the Control Bottles. Add 1 drop of the Negative or Positive Control into the Tube.
- Place a clean Swab into the Tube and follow the instructions for testing the patient Swab.

## LIMITATIONS

The contents of this kit are for use in the **qualitative** detection of Group A Streptococcal antigen from throat swab specimens and culture colonies only. Failure to follow the test procedure and interpretation of test results may adversely affect performance and/or produce invalid results.

The test detects both viable and nonviable Group A Streptococci and may yield a positive result in the absence of living organisms.

Respiratory infections, including pharyngitis, can be caused by Streptococcus from serogroups other than Group A as well as other pathogens. The QuickVue Dipstick Strep A Test will not differentiate asymptomatic carriers of Group A Streptococcus from those exhibiting Streptococcal infection.<sup>6</sup>

Some commercial controls may contain interfering additives and are not recommended for use in the QuickVue test.

Test results must always be evaluated with other data available to the physician. A negative test result might occur if the level of extracted antigen in a sample is below the sensitivity of the test or if a poor quality specimen is obtained. Additional follow-up testing using the culture method is recommended if the QuickVue test result is negative.

## EXPECTED VALUES

It is believed that approximately 19% of all upper respiratory tract infections are caused by Group A Streptococci.<sup>7</sup> Infection is most prevalent in winter and early spring, with most cases arising in patients living in highly populated areas.

## PERFORMANCE CHARACTERISTICS

### Field Study

In a multi-center field evaluation, a total of three hundred twenty-nine (329) throat swab specimens were collected from patients presenting with pharyngitis. Each Swab was inoculated onto a sheep blood agar plate, then tested by the QuickVue Dipstick Strep A Test. Plates were incubated for 24 to 48 hours at 37°C with a Bacitracin disk. Presumptive GAS colonies were confirmed with commercially available Strep A testing kits.

Of the three hundred twenty-nine (329) total specimens, one hundred ninety-two (192) were negative and one hundred thirty-seven (137) were positive by culture. The QuickVue test identified one hundred eighty-eight (188) of the culture negative and one hundred twenty-six (126) of the culture positive for a specificity of 98% and a sensitivity of 92%. The 95% confidence intervals were calculated to be 95% to 99% for specificity and 86% to 96% for sensitivity. Overall agreement between culture and the QuickVue Dipstick Strep A Test was 95%. These study results demonstrate that no statistical differences were observed between the QuickVue test and standard culture techniques.

	Culture +	Culture -	
QuickVue +	126	4	<b>Sensitivity:</b> 126/137 = 92%
QuickVue -	11	188	<b>Specificity:</b> 188/192 = 98%
<b>Total</b>	137	192	<b>Agreement:</b> 314/329 = 95%

Culture Classification	QuickVue Results
Total Negative <b>Specificity</b>	188/192 (98%)
Rare	1/3 (33%)
1+(≤ 10 colonies)	16/20 (80%)
2+(11-50 colonies)	46/49 (94%)
3+(> 50 colonies)	28/30 (93%)
4+ (predominant growth)	35/35 (100%)
Total Positive <b>Sensitivity</b>	126/137 (92%)
<b>Total Agreement</b>	314/329 (95%)

In addition, the QuickVue test was used to confirm the identification of Group A Streptococcus on blood agar plates. As a culture confirmation, the test was 100% sensitive.

## Cross-Reactivity

The following organisms tested at levels of approximately  $1 \times 10^7$  organisms/test and greater were all found to be negative when tested in the QuickVue test.

Streptococcus Group B	<i>Enterococcus faecalis</i>
Streptococcus Group C	<i>E. coli</i>
Streptococcus Group F	<i>Hemophilus influenza</i>
Streptococcus Group G	<i>Klebsiella pneumoniae</i>
<i>Streptococcus pneumoniae</i>	<i>Neisseria gonorrhoea</i>
<i>Streptococcus mutans</i>	<i>Neisseria meningitidis</i>
<i>Streptococcus sanguis</i>	<i>Neisseria sicca</i>
Streptococcus ungrouped	<i>Neisseria subflava</i>
<i>Branhamella catarrhalis</i>	<i>Pseudomonas aeruginosa</i>
<i>Bordetella pertussis</i>	<i>Staphylococcus aureus</i>
<i>Candida albicans</i>	<i>Staphylococcus epidermidis</i>
<i>Corynebacterium diphtheria</i>	<i>Serratia marcescens</i>

## POL STUDIES

An evaluation of the QuickVue test was conducted at three Physicians' Offices using a panel of coded specimens. Testing was performed by physician office personnel with diverse educational backgrounds and work experiences at three different locations. The proficiency panel contained negative, low positive and moderate positive specimens. Each specimen level was tested at each site in replicates of at least five over a period of 3 days. No significant differences were observed within run, between runs or between sites.

## ASSISTANCE

If you have any questions regarding the use of this product or to report a product problem, please contact Quidel Technical Support at 1.800.874.1517 (in the U.S.) or [technicalsupport@quidel.com](mailto:technicalsupport@quidel.com). If outside the U.S., further information can be obtained from your distributor, or directly from Quidel at one of the numbers listed below. Reference [quidel.com](http://quidel.com) to see more options for Support.

Country	Phone	E-Mail Address
Europe, Middle East and Africa	+353 (91) 412 474 (main) 1800 200441 (toll free)	<a href="mailto:emeatechnicalsupport@quidel.com">emeatechnicalsupport@quidel.com</a>
Austria	+43 316 231239	
Belgium	+32 (2) 793 0180	
France	0 (805) 371674	
Germany	+49 (0) 7154 1593912	
Netherlands	0 800 0224198	
Switzerland	0 800 554864	
United Kingdom	0 800 3688248	
Ireland	+353 (91) 412 474	
Italy	+39 (800) 620 549	
North America, Asia-Pacific, Latin America	858.552.1100	<a href="mailto:technicalsupport@quidel.com">technicalsupport@quidel.com</a>
Canada	437.266.1704 (main) 888.415.8764 (toll free)	<a href="mailto:technicalsupport@quidel.com">technicalsupport@quidel.com</a>
China	0400 920 9366 or +86 021 3217 8300	<a href="mailto:chinatechnicalservice@quidel.com">chinatechnicalservice@quidel.com</a>





## REFERENCES

1. Youmans G.P., Paterson P.Y., and Sommer H.M. Upper Respiratory Tract Infection: General Considerations, in the Biological and Clinical Basis of Infectious Disease, W.B. Saunders Co., Philadelphia, 177–183, 1980.
2. Youmans G.P., Paterson P.Y., and Sommers H.M. The Biological and Clinical Basis of Infectious Disease, 1975, pp. 172–185.
3. Recommendations for the Prevention of HIV Transmission in Health Care Settings, Morbidity and Mortality Weekly Report, Centers for Disease Control, August, 1987.
4. Biosafety in Microbiological and Biomedical Laboratories (BMBL), 4<sup>th</sup> Edition, pp. 15-16. U.S. Department of Health and Human Services, COC, NIH, Washington, D.C. (1999). <http://www.cdc.gov/od/ohs/pdffiles/4thBMBL.pdf>
5. Miller JM and Holmes HT. Manual of Clinical Microbiology, Murray P.R., Baron E.J., Tenover F.C., and Tenover F.C., and Tenover F.C., and Tenover F.C. ed. 6<sup>th</sup> Edition, 1995, 3, 28.
6. Rammelkamp C.H., Jr. *Principles of Internal Medicine*, 8th Edition, 1977, pp. 814-820.
7. Lauer B.A., Reller L.B., and Mirrett S. J Clin Microb. 17:338-340, 1983.

**REF**

20108 – QuickVue Dipstick Strep A 50-Test Kit

20108IN – QuickVue Dipstick Strep A 50-Test Kit (Not available for sale in the U.S.)

20125 – QuickVue Dipstick Strep A 25-Test Kit (Not available for sale in the U.S.)

**IVD**



MDSS GmbH  
Schiffgraben 41  
30175 Hannover,  
Germany



**Quidel Corporation**  
10165 McKellar Court  
San Diego, CA 92121 USA  
**quidel.com**

---

### Swab



MDD 93/42/EEC



Emergo Europe  
Prinsessegracht 20, 2514 AP  
The Hague, The Netherlands



Puritan Medical Products Company LLC  
31 School Street  
Guilford, Maine 04443-0149 USA

**1053411EN00 (02/22)**

---

**REF**

Catalogue number



CE mark of conformity

---

**EC REP**

Authorized Representative  
in the European Community

**LOT**

Batch code

---



Use by



Manufacturer

---



Temperature limitation



Intended use

---

**R<sub>x</sub> ONLY**

Prescription use only



Consult instructions for use

---

**IVD**

For *In Vitro* diagnostic use



Contains sufficient for XX determinations

---

**CONT**

Contents/Contains

---