



Sofia²

Flu+SARS Antigen FIA

Flu A, Flu B + SARS-CoV-2 at Your Fingertips

One Test, Three Results, 15 Minutes



Accurate, objective and automated results in **15 minutes**



Integrated data management automatically stores test and user history



Excellent performance compared to viral culture or molecular methods



User and patient ID captured with **onboard barcode scanner**




Flexible, dual mode testing for high throughput in variety of laboratory environments




Sofia 2 Flu + SARS Antigen FIA - Procedures

Nasal/Nasopharyngeal Swab Procedure

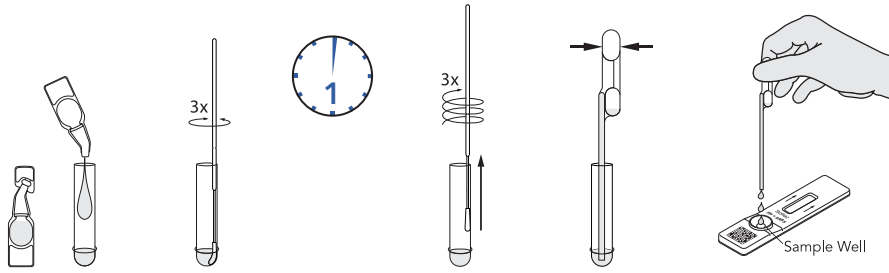


Sofia 2 Development Modes

WALK AWAY Mode – Walk away and multitask

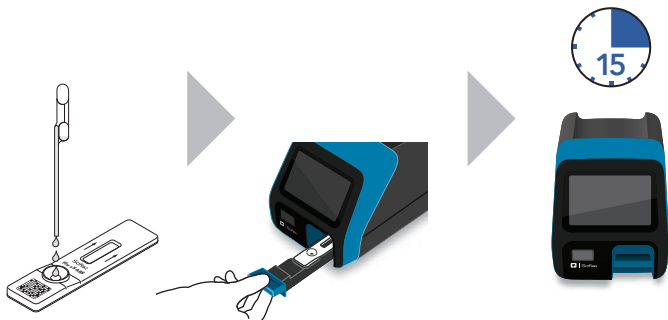

READ NOW Mode – Batch multiple samples per hour

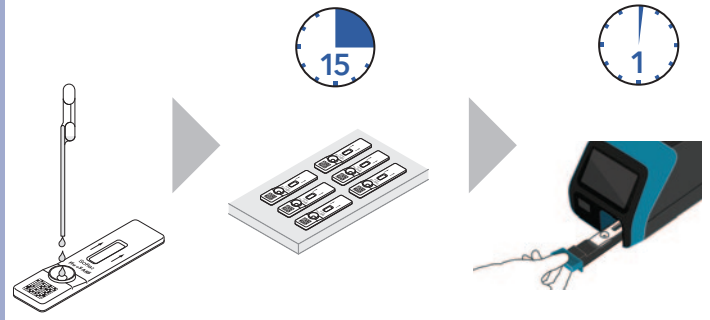

Sofia 2 Flu + SARS Antigen FIA*

- Nasal and nasopharyngeal specimens
- Accurate detection with direct samples
- Results in 15 minutes
- Kit includes 25 sterile nasal or nasopharyngeal swabs and a set of positive and negative control swabs
- Room temperature storage
- 24-month shelf life from date of manufacture

Sofia 2 Flu + SARS Antigen FIA Clinical Performance*

	Sample Type	Sens/PPA*	Spec/NPA*
Influenza A	Nasal	90%	95%
	Nasopharyngeal	97.1%	94.6%
Influenza B	Nasal	89%	96%
	Nasopharyngeal	90%	97%
SARS-CoV-2*	Nasal	95.2%*	100%*

*Results for SARS-CoV-2 are expressed as PPA/NPA and were tested on nasal swabs in the study, however, may also be used with NP swabs.

Refer to the Package Insert for additional performance claims.



Sofia 2 Flu + SARS Antigen FIA – 25 Test Kit: Catalog #20380

SS2038002EN00 (09/21)



EMEA HQ: Quidel Ireland Ltd.
2nd Floor, Merchants Square
Merchants Road Galway,
Ireland H91 ETN2



+44 800 368 8248



1800 200 441

EMEA: +353 (91 412 474

emeacustomerservice@quidel.com